




# Foam Tipped Applicator

Acetal handle

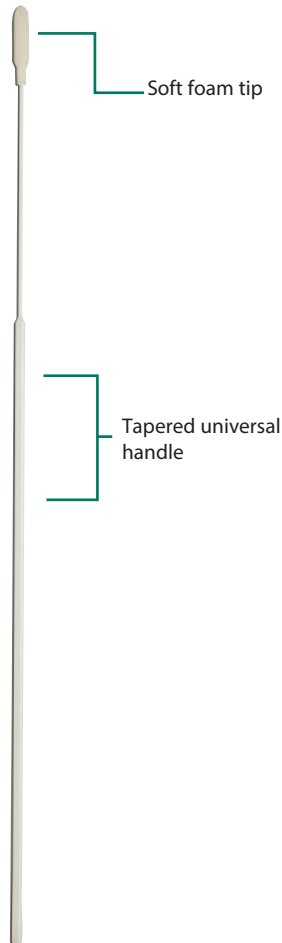

## Description

This specialty applicator is tipped with 100ppi medical grade polyurethane foam. It is bonded securely to a tapered acetal handle using a thermoplastic adhesive and heat stake process. Often used as a general swab for medical or veterinary diagnostics.

## Dimensions

6" Overall Length (152.4 mm)

### Tip

- 0.141" Dia. (3.581 mm)
- 0.547" Length (13.894 mm)

### Handle

- 0.098" Dia. (2.489 mm) tapered to 0.049" (1.244 mm)
- 5.996" Length (152.298 mm)

## Qualities

- Medical grade quality
- Foam has excellent collection and release properties
- Soft foam tip combined with flexible handle provides maximum patient comfort
- Made in USA

## Applications

- Specimen collection
- Nasopharyngeal or Ureathral specimen collection
- Recovery or micro-organisms and fluid
- General swab for medical or veterinary diagnostics
- Bacterial and virology

Ref #	Packaging
25-1406 1PF 50 Sterile	10/50/1 paper/poly pouch (10 boxes of 50 individually wrapped packages per case, 1 swab per/pkg)



**Puritan®**