




Polyester Tipped Applicator

Tapered mini tip, Wood handle

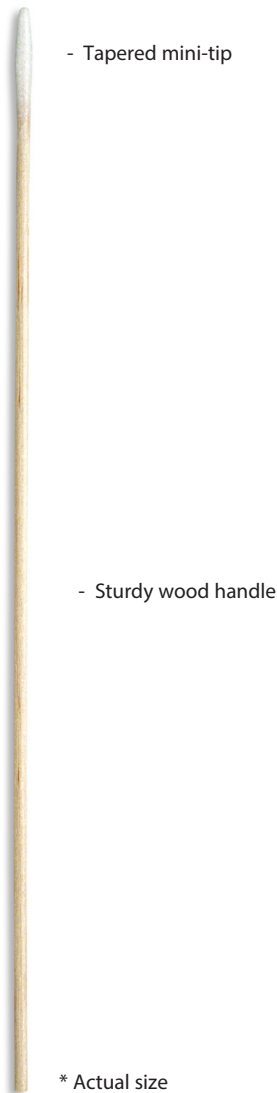

- Tapered mini-tip

- Sturdy wood handle

* Actual size

Description

This synthetic tipped applicator is made of medical-grade polyester fiber that is carded in our factory to maintain consistency. It contains no fatty acids or latent sugars, making it ideal for specimen transport, lab work, and diagnostic kits. It is bonded securely to a sturdy wood handle by an aqueous based adhesive.

Dimensions

6" Overall Length (130.175 mm)

Tip

- 0.12" Dia. (3.048 mm)
- 0.5" Length (12.7 mm)

Handle

- 0.083" Dia. (2.108 mm)
- 5.812" Length (147.625 mm)

Qualities

- Medical grade quality
- Polyester tip contains no fatty acids or latent sugars
- Tapered mini-tip for hard-to-reach or small areas
- Sturdy wood handle, naturally biodegradable
- Made in USA

Applications

- Specimen collection
- Patient care
- Surface collection

Ref #	Packaging
25-826 2WD Sterile	10/100/2 (10 boxes per case of 100 individually packages per box, 2 p/pkg)



Puritan®